


# WOW

## Carrying sheet with 10 handles

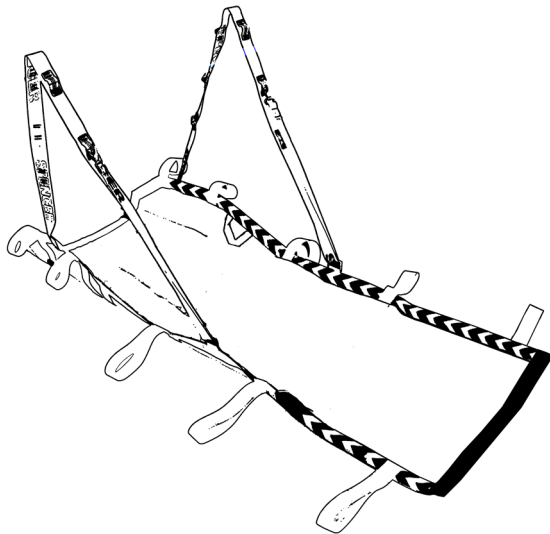

The carrying sheet WOW, is a device intended to transfer patients from the stretcher to the bed or for transport from the rescue place to the stretcher for patients who do not have spinal injuries.

### Specific Characteristics

- Manufactured with fireproof materials
- Easy to clean thanks to the use of non-absorbent fabrics
- Reinforced nylon webbing to ensure high resistance and durability
- The width of the tape used for the handles, ensures an easy grabbing increasing safety
- On the patient side there are handles dedicated to the patient that make him more collaborative and increase safety
- The special shape, offers high lateral containment during transport
- The supplied adjustable shoulder straps ensure a better distribution of the load on the rescuer
- The longitudinal aluminium rods, are removable through the dedicated holes and can be folded into 3 parts. When inserted, they increase the rigidity of the sheet, allowing to distribute the load over a wider surface and allowing to limit the flexion of the device
- The device allows also the transport of the patient in sitting position
- Requires low storage space

### Compliant with

UNI EN 1865-1

### Technical data

Length	1870 ± 10 mm
Upper part width	960 ± 5 mm
Lower part width	700 ± 5 mm
Storage dimensions <sup>(1)</sup>	500x600x150 mm
Dimensions of a folded rod	75x450x30 ± 5 mm
N° of handles	10
Width of tapes	50 mm
Upper handles length	160 mm
Lower handles length	260 mm
Rod length	1210 ± 10 mm
Loading capacity	150 kg
Materials	Nylon, PVC, Al
Weight with rods	4,2 ± 0,2 kg
N° of shoulder straps adjustments	6

<sup>(1)</sup> The dimension is referred to the wrapped sheet. It is indicative, since the device can be folded in many ways.

Class I Med.Device according to Dir. 93/42/CEE

WOW - CARRYING SHEET 10

HANDLES

QM40601B

CND Classification



V0899

Rev.1 (31/03/2016)

UNCHECKED COPY – further revisions will be available on <http://support.spencer.it>

Spencer Italia S.r.l. a socio unico Collecchio (PR) Italia Tel. +39.0521.541111 Fax +39.0521.541222

[www.spencer.it](http://www.spencer.it)