


WOW

Carrying sheet with 10 handles

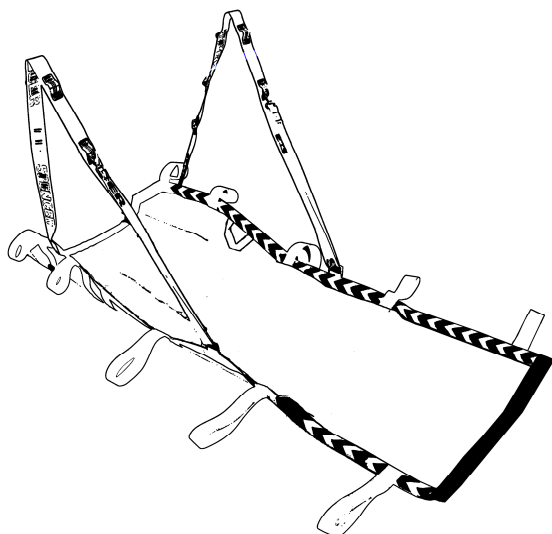

The carrying sheet WOW, is a device intended to transfer patients from the stretcher to the bed or for transport from the rescue place to the stretcher for patients who do not have spinal injuries.

Specific Characteristics

- Manufactured with fireproof materials
- Easy to clean thanks to the use of non-absorbent fabrics
- Reinforced nylon webbing to ensure high resistance and durability
- The width of the tape used for the handles, ensures an easy grabbing increasing safety
- On the patient side there are handles dedicated to the patient that make him more collaborative and increase safety
- The special shape, offers high lateral containment during transport
- The supplied adjustable shoulder straps ensure a better distribution of the load on the rescuer
- The longitudinal aluminium rods, are removable through the dedicated holes and can be folded into 3 parts. When inserted, they increase the rigidity of the sheet, allowing to distribute the load over a wider surface and allowing to limit the flexion of the device
- The device allows also the transport of the patient in sitting position
- Requires low storage space

Compliant with

UNI EN 1865-1

Technical data

Length	1870 ± 10 mm
Upper part width	960 ± 5 mm
Lower part width	700 ± 5 mm
Storage dimensions ⁽¹⁾	500x600x150 mm
Dimensions of a folded rod	75x450x30 ± 5 mm
N° of handles	10
Width of tapes	50 mm
Upper handles length	160 mm
Lower handles length	260 mm
Rod length	1210 ± 10 mm
Loading capacity	150 kg
Materials	Nylon, PVC, Al
Weight with rods	4,2 ± 0,2 kg
N° of shoulder straps adjustments	6

⁽¹⁾ The dimension is referred to the wrapped sheet. It is indicative, since the device can be folded in many ways.

Class I Med.Device according to Dir. 93/42/CEE

WOW - CARRYING SHEET 10
HANDLES

QM40601B

CND Classification



V0899

Rev.2 (18/07/2017)

UNCHECKED COPY – further revisions will be available on <http://support.spencer.it>
Spencer Italia S.r.l. Sala Baganza (PR) Italia Tel. +39.0521.541111 Fax +39.0521.541222

www.spencer.it