


## BOB MASTER

### Stretcher support

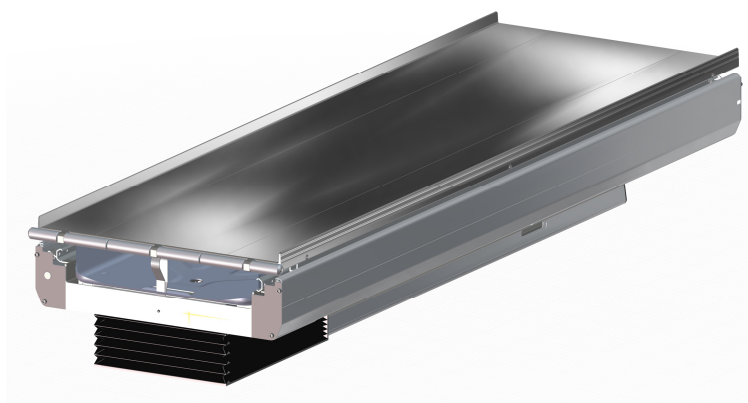

Stretcher support intended for installation on ambulances. The device is able to assume several positions that facilitate both loading and unloading of the stretcher, allowing at the same time easier access to the stretcher.

#### Specific features

The structure of the support plane makes it non-deformable  
 Loading system with roller. Thanks to the low height that the device can reach, a slide is not needed  
 Dedicated compartment to store immobilization systems like spinal boards or pick-up stretchers.  
 Can assume the positions tilt front and tilt back  
 Adjustable loading height to make the stretcher loading easier  
 Extruded aluminum board easy to clean  
 The control panel is integrated in the structure. Additional control panel can be installed, for example on the ambulance wall  
 Operating handles are easily accessible  
 Suspension automatically regulated according to the patient weight  
 Lockable allowing to perform CPR  
 Extraction in tilted position to allow easy loading/unloading of the stretcher  
 Can be moved sideways manually (both left and right), extracted and moved forward  
 The device automatically returns in horizontal position once it is moved inside the vehicle

#### Reference standard

10g tested according to EN1789:2007+A2:2014  
 EN 1865-5:2012  
 ECE-R-10

#### Optional accessories

BO10107 – Basin for spinal board housing  
 EXTERNAL CONTROL PANEL

#### Technical data

Length	1990 mm
Width	600 mm
Minimum height	290 mm
Maximum height	390 mm
Sideways movement	230 ± 10 mm
Extraction tilted (approx)	1120 mm
Tilt front angle	Max 4°
Tilt back angle	Max 15°
Loading height <small>referred to the installation floor</small>	Min (-9cm, max +29cm)
Materials	Al, Steel, PE, Nylon
Loading capacity	300 kg
Weight	125 ± 2 kg

#### Device in compliance with dir. 2006/42/CE

BO10008C – BOB MASTER SHOCK  
 ABSORBING STRETCHER SUPPORT  
 SIDEWAYS MOVABLE AND  
 EXTRACTABLE



Rev.0 (27/09/2018)

UNCHECKED COPY – further revisions will be available on <http://support.spencer.it>  
 Spencer Italia S.r.l. Sala Baganza (PR) Italia Tel. +39.0521.541111 Fax +39.0521.541222  
[spencer.it](http://spencer.it)